



DIRECT HITS

Keys to Successful Selling

North American Industry Classification System

325412 – Pharmaceutical Preparation Manufacturing
541710 – Biotechnology Research and Development

Application:

Disposable tube sets used to transfer cell culture media between seed train bioreactors.

In This Application

The SaniQuik™ Connection provides a quick and easy link between stainless equipment and disposable tube sets. SaniQuik's integral sanitary termination attaches directly to the bioreactor with a sanitary gaskets and tri-clover clamp. Once attached to bioreactor the connection's integral coupling adaptor provides the interface to Colder's MPC and MPX coupling used on the disposable tube set.

Media:

Cell Culture

Pressure:

Up to 60 psi

Temperature:

Ambient


CPC Part Numbers:

SQCC221212M, SQCC222424M
SQCX221212M, SQCX222416M, SQCX222424M

Features/Benefits:

- 1) Provides a quick and easy link between disposable systems and existing stainless steel processing equipment
- 2) Connects to hard-plumbed systems with standard sanitary gasket and tri-clover clamp
- 3) Mates with Colder's MPC and MPX disposable couplings

Tube Sets

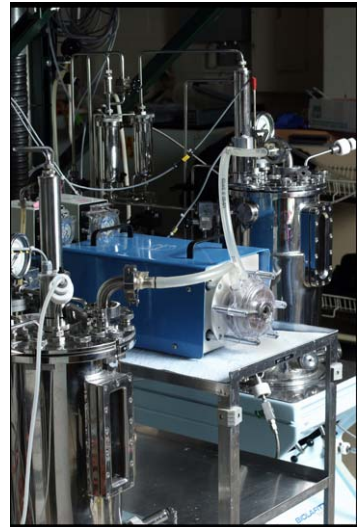

Features/Benefits:

- 4) Increases operational efficiency
- 5) Withstands multiple sterilization cycles
- 6) Ideal for use in feeding, cell harvesting, and sampling applications

All sales are subject to Colder Products Company's limited express warranty set forth in the Colder catalog. Contact your local distributor or Colder Customer Service for warranty provisions.

1001 Westgate Drive, St. Paul, MN 55114 U.S.A. • Phone 651-645-0091 • Fax 651-645-5404 • www.colder.com

NAICS
325412
541710

Direct Hit
BIO-2